

World recognizes the technology and effects!

Lucid Korea's Ortho-K lenses' (LK -Lens, Toric Premier Lens) technology and its effects have been recognized worldwide and it has been sold to many countries abroad such as China, India, Ukraine, Vietnam and so forth.

Approved by KFDA

MDR CE approval is expected by end of this year.

LK -Lens and Toric Premier Lens have been **approved by KFDA**, and in Korea, we are the first one who got approved.



Effects on controlling myopia progress

According to international ophthalmology journals IOVS (April, 2011), PLOS (April, 2015), there was meaningful change on axial length compare to the comparative group, it shows effects on controlling myopia progress.



Tel +82-2-2282-6011

B-1401, Hangang XI tower, 401 Yangcheon-ro, Gangseo-gu, Seoul, Korea

E-mail claire@lucid.co.kr, hyein@lucid.co.kr

Ortho-K

LK - Lens Toric Premier Lens

Ortho-K Lenses are to correct your vision while sleeping,
it can make you see clearly without lenses or glasses all day long

www.lucid.co.kr




Lucid Korea's Ortho-K Lenses

Get liberated from glasses and lenses!


Ortho-K Lenses, how can you see well without wearing lenses during the daytime?

Lucid Korea's Ortho-K Lenses' specialized structure enables you to see well without lenses and glasses during daytime by pressing cornea while sleeping. If you stop wearing lenses, then your cornea goes back to the original shape, so it doesn't prevent you from getting a Lasik surgery or damage your eye health. Ortho-K Lenses are known as 'Dream Lenses', cause it can only require you to wear them for about 6 hours of your sleep and you can see well during the daytime.

Before wearing Ortho-K Lenses

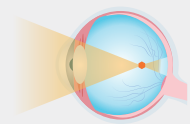


After wearing Ortho-K Lenses

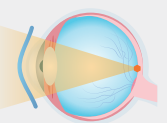

WHY LUCID KOREA?

5 curves structure with better safety and better effects

Compared to 4 curves, Lucid Korea's Ortho-K lenses have 5 curves which minimize the lens-to-cornea contact, it enhances patient comfort, safety and accuracy of vision correction.

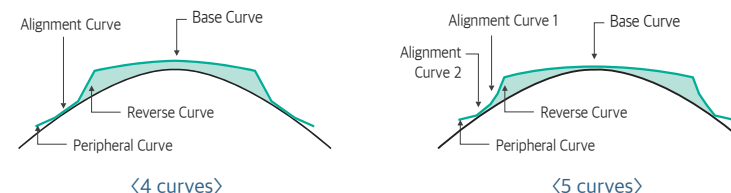

30 years of research and expertise in Ortho-K

Since its establishment in 1996, Lucid Korea has spent the past 30 years researching, developing, and refining Ortho-K lens technology.

Accurate and safe fitting procedure, good vision correction effect

With just one fitting procedure, it can correct vision for myopia and astigmatism, even with severe myopia up to -6.00D, it is very effective.

Ortho-K LK Lenses, What part is good?



SAFE



Good for children and teenagers!

For children who feel uncomfortable wearing lenses or glasses, or afraid of getting Lasik surgery, Lucid Korea's Ortho-K lenses are great choice for them. It enables them to see well during the daytime without glasses or lenses!

No age limit!

Regardless of age and gender, anyone can try Lucid Korea's Ortho-K lenses. For children, parents' guidance is required.

Side effects free!

When your eyes feel uncomfortable, you can simply stop wearing lenses, and the cornea goes back to original shape. So it doesn't have side effects.


Ortho-K LK, TORIC Lens Specifications

- Myopia up to -6.00D, and Astigmatism up to -4.00D
- Dk = 100
- Flat K 7.50~8.39
- Diameter 10.0~11.2mm (0.2mm increments)


Made With



Lucid Korea's Ortho-K lenses are manufactured with Bausch & Lomb's 'Boston XO' which has high oxygen permeability (Dk = 100).

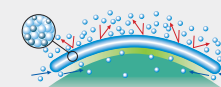

PMMA LENS
(Dk=0)

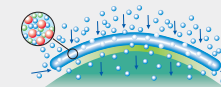

LUCID KOREA'S LENS
(Dk=100)